

medITECH

MEDICAL GAS EQUIPMENT

Twin Flow Selector Wall Mounted

With 2M supply hose **Also available with Rail Mount**

- No Calibration required EVER
- Lightweight construction
- Available with International Fittings
- 7 Year Warranty • 7 Year service Interval




Excellent performance and robustness, coupled with a very lightweight design. For use in all Hospital and Emergency environments, the Twin Flow Selector is designed to be plugged into Medical Oxygen Wall Outlets, whilst being attached directly to the wall, or attached to a rail. The Meditech Twin Flow Selector has been designed to deliver precision flow rates, without the need for any calibration. It's fitted with either a universal wall mounting plate or a universal size rail mount for most rail systems. It has a standard 2M hose length.
(other hose lengths available upon request).


Can be supplied in either of two formats:

- Push-Fit 'Firtree' Fitting for direct attachment on an Oxygen Therapy delivery system
- With a threaded outlet fitting for the attachment of a humidifier.
(Humidifiers to be purchased separately)

PROBE TYPES:

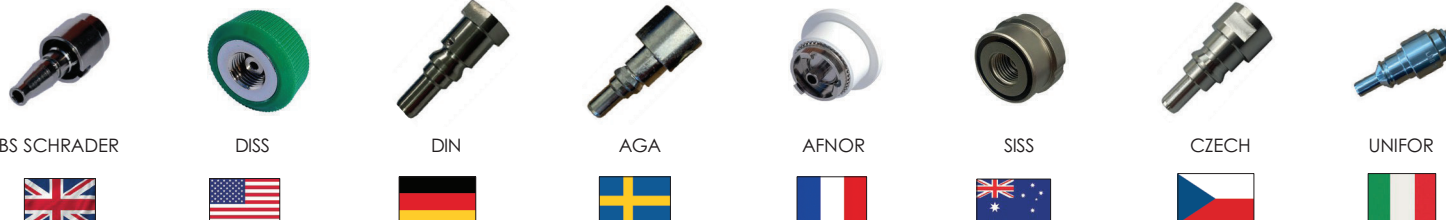

THERAPY TYPES:

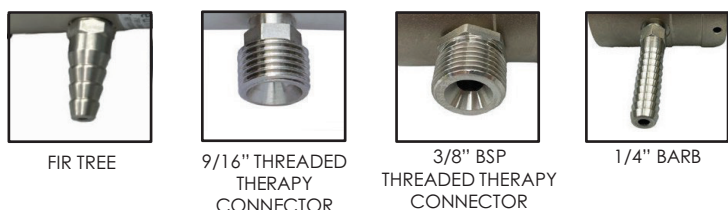

AVAILABLE GAS TYPES:

- O₂
- MA4

SPECIFICATION

Intended Gas:	Medical Oxygen (O ₂ , MA4)
Inlet pressure:	4 bar (400 kPa)
Set Flow Rates in lpm:	Three flow rates options available: 1/2, 1, 2, 3, 4, 6, 8, 10, 15 lpm 1/2, 1, 2, 3, 4, 6, 8, 10, 15, 25 lpm 10, 20, 30, 35, 40, 50, 55, 60, 70 lpm
Inlet Connections:	Schrader (BS 5682:2015) and other international connections
Outlet (Patient) Connections.	Firtee style barbed outlet, 9/16 threaded outlet, 3/8" BSP threaded outlet & 1/4" barb
Intended Life:	10 Years
Service Interval:	7 Years
Warranty Period:	7 Years
Complies With:	ISO 15002: 2023

TO PURCHASE OR FOR FURTHER DETAILS:
Contact us by email or telephone for more details,
including specific part numbers, prices, etc.

B.N.O.S. Meditech Ltd.
9 Fifth Avenue, Bluebridge Industrial Estate,
Halstead, Essex CO9 2SZ, England, UK.
Tel: +44 (0) 1787 479475 • Fax: +44 (0) 1787 477747
E-mail: sales@meditech.uk.com • www.meditech.uk.com



Meditech Flow Selectors are CE marked in accordance with the Medical Devices Directive in the EU and UKCA to the Medical Devices Regulations 2002. All items in this range are sold with a return to manufacturer warranty against manufacturing defects. These products are part of a series of medical gas and resuscitation products available from B.N.O.S. Meditech. If the product you require is not shown here, or if you require more information, please contact your distributor or contact us directly. Specifications may be subject to change without notice. B.N.O.S. Meditech Ltd Terms and Conditions of sale apply, available on application or from our website. Entonox® is a trade mark of the Linde Group

medITECH
MEDICAL GAS EQUIPMENT

Distributor:

medITECH
MARINE

E: sales@meditechmarine.com

Doc 778M Issue D 18/05/26